OPERATING TABLE ACCESSORIES TAB Series





MANUFACTURING. DEVELOPMENT.CUSTOMIZATION

Qualität ist unsere Priorität

International group medical technology and electronics GmbH

Anesthesia Frame & Binding Strap & IV Rod



L type Anesthesia Frame

Use to place curtain It's composed of stainless steel anesthesia frame and clamp welded together

Main features

Use high quality su304 stainless steel

Fixed type, no need equip with clamp



Body Binding Strap

Usetofixpatient'sbody It'scomposedofbinding strap and fixed

Main features

lse high quality su304 stainless steel

Size: 520x160 (mm);

Length of Strap: 1140mm

Install on the side rail



Standard Anesthesia Frame

- · Use to place curtain
- Able to adjust the height, The curtain frame is stretchable

Main features

Use high quality su304 stainless steel

Install on the side rail, clamp sold separatly



Arm hanging belt

Usetohangpatient's arm

It's composed of hanging belt and buckle It can be adjusted up-down, left-right

Main features

Use high quality su304 stainless steel
Install on the side rail, clamp sold separatly



Luxurious Anesthesia Frame

- It's able to make support surface through supporting arm
- ' The height is adjustable

Main features

Use high quality su304 stainless steel

Install on the side rail, clamp sold separatly



IG Tab 21 -1601.0150

Hand Binding Strap

Use tofixpatient'sarm

It's composed of hand binding strap and support buckle

Main features

se high quality su304 stainless stee

Install on the side rail



Adjustable IV Rod

- · Swingable IV rod
- · Use to hang infusion bags and bottles
- · It can swing freely, height is adjustable

Main features

Use high quality su304 stainless stee

Adjustable height range 800mm-1500mm

Install on side rail, including clamp



IG Tab 21 -FJ-36

Leg Binding Strap

Use to fixpatient'sleg

It's composed ofbinding straps for leg and leg board

Main features

se high quality su304 stainless stee

nstall on the side rail



Hand Support & Leg support



Arm Rest (Low Position)

- · Usetosupportandfixpatient's arm
- It'scomposedofarmrest,arm rest adopter,adjustableclampand pad
- * Theangelcanbeadjustedby itself
- Itcanbeadjustedthewidth and height vertically and horizontally

Withtwohandbindingtraps to fix arm

Main features

Ise high quality su304 stainless steel

Inside of pad is soft sponge

Thickness of pad: 60mm

Install on side rail, including two hand binding straps, clamp sold separatly



10 1ab 21 -1310.1000

- Extension Board (Pelvis)

 Mainlyusedtoextendtheelectric operating table
- It'scomposedofRadiolucent board and supportframe

Main features

o use high quality su304 stainless steel and

Apply to COMFORT 300 operating table

Extension Board Dimension: 524x250 (mm)



Arm Rest (Lateral Position)

- $\bullet \ \ {\sf Useto support and fixpatient's arm}$
- It'scomposedofsupportingrod,
- supporting board and pad
- Theangelcanbeadjusteditself, and the height as well.

Withtwohandbindingtrapsto fix arm

Main features

lee high quality su30/, stainless steel

nside of pad is soft sponge

Thickness of pad: 60mm

Install on side rail, including two hand binding straps, clamp sold separatly



Hand Protector

- Mainlyusetoprotect patient's body and arm
- ' Usehighqualitystainless steel

Main features

side of liner is soft sponge

Dimension: 460x150x130 (mm)

Install on the side rail



Height Adjustable Arm Rest

- $\bullet \ \ {\sf Useto support and fixpatient's arm}$
- It'scomposedofsupportingrod,
- supportingboardandpad
- Itcanrotatefreely,heightisadjustable
 Withtwohandbindingtrapstofix arm

Main features y su 304 stainless stee

nside of pad is soft sponge

Thickness of pad : 60mm

Lxw: 660x140 (mm)

Install on side rail, including two hand binding straps, clamp sold separatly



Luxurious Lithotomic Leg Support

- · bigger size
- · adjustable freely
- · with binding strap

Main features

Use high quality su304 stainless steel Its

nonge with leather soft and cozy Install

on side rail, clamp sold separatly



Arm Surgery Assistant Board /Carbon Fiber Surgery Board

- ·Use to support and fix patient's arm
- · It's composed of supporting frame, supporting board and pad
- Radiolucent

Main features y su304 stainless steel and radiolucent board (Carbon fiber is optional)

Inside of pad is soft spond

Thickness of pad: 60mm

LxW: 800x350(mm)

Install on side rail ,need equip with another tow clamps(#1402.2000)



Leg Hanging Frame

- Usetosupportandfixpatient's leg andhand
- it'scomposedofsupporting frame and hangingstrap
- Heightisadjustable
- · Itrotatefreely

Main features

Ise high quality \$1304 stainless steel

nstall on side rail, clamp sold separatly



Skull Clamp & Head Support & Leg Plate



Head Support for Ophthalmology Surgery

- · Usetosupportpatient'shead
- meanwhilesupportsurgeon'sarm
- It's composed of support frame, head support and radiolucent board
- Theheightofsupportframeisadjustable

Main features

high quality su304 stainless steel

Inside of head support is soft sponge



Head Support Adaptor

- · Usetoequipwith1514.0000head support
- It'scomposedofsupport frame and lockinghandle
- Providedifferenttypeof head support interchangeonoperating table

Main features

Use high quality su304 stainless steel

1413.1200 (for YS-3)

1413.0100(for 300)



Skull Clamp (Neurosurgery) ·

- Usetofixpatient'sheadwhenoperate
- neurosurgery
- It'scomposedofMainbody,skull fix
- frameandspurshapeheadsupport Able to rotate freely

Main features

se high quality Aluminum

Install on the bridge support (Neurosurgery support)



Cylinder Buttock Support ·

 Usetosupportpatient'sbuttock and leg It'scomposedofsupportframe and buttock support

Main features

se high quality su304 stainless steel

Inside of buttock support is soft sponge

Install on the side rail, clamp sold separatly



Bridge Support (Neurosurgery support)

- It support skull clamp to fix patient's
- ·It's composed of bridge support and locking handle
- •able to adjust the height

Main features

Jse high quality su304 stainless steel

Install on the side rail, clamp sold separatly



Leg plate

Usetosupport patient's leg
It'scomposed of support frame and leg

Main features

Jse high quality su304 stainless steel

Inside of leg plate is soft sponge

Liner pad dimension: 280x270x60 (mm) Install on the side rail, clamp sold separatly



Double Joint Head Support

- · Usetosupportpatient's head
- It'scomposedofsupportingframe and
- headsupport

Abletofreelyadjusttheheightof joint

Main features

Ise high quality su304 stainless steel

Inside of head support is soft sponge

its be used with 1413.0100/1413.1200 head support adopter



Leg Plate (Pair)

Usetosupportand fix patient's leg It'scomposedofsupporting frame and linerpad

abletorotatefreely

Main features

Use high quality su304 stainless steel

Inside of liner pad is soft sponge

Liner pad dimension: 280x270x60 (mm)
Install on the side rail, including clamp



Body Support & Shoulder Support & Accessories Cart



IG Tab 21 -1601.0125

Economic body support

- · · Usetosupportorfixpatient's body
- . It's composed of body support and

The height can be adjusted

Main features

High quality su304 stainless steel Inside of liner is soft sponge Liner size : 215x100x60(mm) Install on the side rail, clamp sold separatly



Standard shoulder & body support

- ·Use to support or fix patient's shoulder
- •The height can be adjusted •360° rotation is adjustable

Main features

High quality su304 stainless steel Inside of liner is soft sponge Liner size : 215x100x60(mm) Install on the side rail, clamp sold separatly



IG Tab 21 -1601.0315

Luxurious shoulder support

- $\cdot \cdot$ Use to support or fix patient's shoulder
- It'scomposedofshouldersupport and support rod
- The height can be adjusted

Main features

clamp sold separatly

High quality su304 stainless steel Inside of liner is soft sponge Liner size : 215x100x60(mm) Install on the side rail,



Shoulder Support

- Use tosupport or fixpatient's shoulder and body
- •The height can be adjusted

Main features

igh quality su304 stainless steel

Easy to clean disinfect

Length300mm,Width 300mm,Height 90mm, Load capacity:10Kg



IG Tab 21 -1601.0096-A

Luxurious body support A

- Usetosupportorfixpatient'sshoulder and body
- It'scomposedoflinersupportand support rod
- * The height canbe adjusted
- * 360°rotationis adjustable

Main features

High quality su304 stainless steel Inside of liner is soft sponge Liner size: 215x100x60(mm) Install on the side rail, clamp sold separatly



X-ray Cassette Tray

It'scomposedofX-raycassette and
handle

Main features

Ico cu704 staipless steel

X-Ray Cassette Dimension: 450x400 (mm)



IG Tab 21 -1601.0096-B

Luxurious body support B

- · · Use to support or fix patient's body
- It's composed of body support and
- support rod
- * The height can be adjusted 360° rotation is

Main features

High quality su304 stainless steel

Inside of liner is soft sponge

Liner size: 215x100x60(mm)

Install on the side rail, clamp sold separatly



Accessories Cart

- ·Mainlyuseto placetable accessories
- •It's composed of main frame and mobile
- ·Able to 360° rotation freely

Main features

Use su304 stainless steel

With self locking castors



Clamp



Non-adjustable Clamp

- · Usetofixroundorsquaresurgical
- * It'scomposedofClampand handle
- · Accessoriespositioncannot be adjusted

Main features

Applicable to Φ ≤ 16mm round rod accessorie

Install on the side rail



Adjustable Clamp

Usetofixoperatingtable accessories It'scomposedofClamp and handle It can be adjusted freely

Main features

It can fix Φ≤16mm round rod accessories

Install on the side rail



Square Clamp

Usetofixoperating table accessories It'scomposedofClamp and handle It cannot beadjusted

Main features

Applicable to $\Phi \leq 16$ mm round rod accessories and profile ≤ 6 * 25mm rectangular accessories



IG Tab 21 -1601.0026

Square Pipe Clamp

Usetofixoperatingtable accessories It's composed of clamp and handle

Itcanbeadjustedfreely

Main features

Fix 25mmx25mm square pipe

Install on the side rail

Waste Basin & Sponge Position Pad



Stainless Steel Waste

Use to contain surgery waste

IG Tab 21 -1404.3000



Semi-round Pad (Sponge)

To be used under the knee



Head Ring



Knee Position Pad (Sponge)





Head and Shoulder Pad (Sponge)

IG Tab 21 -1415.0005



Single Side Instrument

Use to place instrument It is used with instrument

basin 1409.0001



IG Tab 21 -1415.0006

Lateral Position Leg

Tunnel Pad



(Sponge)



(Sponge)



Accessories for Standard Configuration & Orthopedic tractor



Arm Rest Hand Binding Strap

 Adjustable byjoint adapotor, flexible and meet different patient's surgical position

Easy to clean and disinfect

Main features Use high quality SU304 stainless steel Thickness of pad:60mm Load capacity:20Kgs



Economic Lithotomic Leg Support

- · Use to support patient's leg
- · It's composed of supporting rod, leg pad
- · it rotate freely

Main features Use high quality su304 stainless steel The leg pad is made of polyurethane foam Install on side rail, clamp sold separatly



Body Binding Strap

• •Use tofix patient'sbody
It's composedofbinding strap
and fixed buckle

Main features Install on the side rail Use high quality SU304 stainless steel



IV Rod

- Itissuitable for hanging infusion bagsand bottles
- · Itcanswing freely

Main features Use high quality SU304 stainless stee Easy to clean and disinfect Load capacity:50Kg

Orthopedic tractor #1102.0000

Main features

The table top elevation displacement of the orthopedic tractor is 350 mm. It can be adjusted according to the elevation displacement of the operating table, which is favorable for the medical personnel to operate and use.

The pelvis support plate and pad and the perineum pillar are both made of nonmetallic material, which will not influence the X-ray shooting. During operation, the pad will not be moved. At the same time, the pelvis support plate and pad can be dismounted in order to radiograph, fix film carrier, clean and to disinfect

The orthopedic tractor can be stretched from 0 to 700 mm.

The maximum carrying capacity is 160 kgs.

The maximum traction displacement is 180 mm and the maximum traction power is 400 N.

The orthopedic tractor can be folded outside and inside in order to meet the needs of all surgeries.

The safety design of the orthopedic tractor is in accordance with the general require for safety of Medical Electric Equipment GB9706.1-2007.

The orthopedic tractor can be used more than 5 years.





Other Accessories



Pneumatic stirrup (pair)


- One hand control. Move freely and easily set at required position. With mattress and bidding strap.
- ·Clamp on the side rail.

Main features

Jse high quality su304 stainless steel

Install on the side rail, clamp sold separatly

Inside of leg support is soft sponge

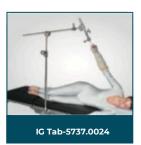

Lateral Shoulder Holder

- · Easy to operate, by simply adjust several traction knob to complete traction setting
- · Simple structure, weightless
- . All components are autoclavable
- · Traction settings from 0-11KGS

Main features

Use high quality su304 stainless steel
Install on the side rail, clamp sold separatly



Economic Precise Leg Support

- · Integrated booth, combine with soft liner
- Able manually adjust and fix to desired position
- · Working load 158KG(350LB)

Main features

th quality su304 stainless steel

Inside of leg support is soft sponge

Install on the side rail, clamp sold separatly



Elbow Surgery Positioner

- · All components are autoclavable
- Quick release allows for intraoperative check for range of motion
- Combines ease of use and adaptability to provide intraoperative positioning for supine elbow surgery

Main features

quality su304 stainless steel

Install on the side rail, clamp sold separatly



Shoulder Surgery Positioner

- Able to quickly and precisely locked patient's arm in desired position
 Single lever releases all joints for positioning and relocks
- Less man needed, simple and easy
- traction capacity: 227KG (500LB)

Main features

se high quality su304 stainless stee

Install on the side rail



IG Tab-5737.0037

Accessories Cart KneeSupport

Mainly use to place table accessories Provide more safety support during surgery It's composed of main frame and mobile Using a disposable conductive non-latex castors

ARabdleiotlou 360° rotation freely

Maaiinn ffeeaattuurreess

પાકેsee csuar3b0o4n s fitablenrl emssa tseterie. With self locking castors

Install on the side rail, clamp sold separatly



Legholder

Braces can safely hold patient's leg 20mm install pole 360° rotation is adjustable

Main features

The Braces is made from soft sponge

Equip with clamp

Install on side rail



Spine Frame

Spine frame provides stable and convenient platform for spine surgeries such as laminectomies, decompressions You can adjust it to the desired lordosis to get optimal surgical site exposure Weight capacity 227KG (500LB)

C arm compatible, use high quality sponge, safe and comfortable

Main features

Using carbon fiber material

Using high quality sponge, safe and comfortable





Knee Positioner

Use to support and fix patient's leg Track Locking Bracket securely locks the base plateto thetable with two thumb

screws Smooth sliding for leg positioning During knees surgery, leg can extent, flex, tilt and rotate controlled by ball joint

Main features

Using aluminum-magnesium alloy material

Disposable sterile boots liner

Install on side rail



Shoulder Surgery Beach Chair

Separate detachable shouldermodule for easy access to shoulder area

The head mounting module can be rotated and tilted

The pneumatic rod can be folded or lowered and fix the control backboard

The two side boards can be detached or part detached

Main features

nside of liner is soft sponge

Install on side rail



COMFORT shoulder support

- · Use to support patient's shoulder
- Ergonomic,specialmodelsuitable for shoulder support
- ' It can be adjusted left and right
- Shoulder support can rotate freely

Main features

High quality su304 stainless steel or aluminun

Inside of shoulder support is soft sponge

Install on the side rail, clamp sold separatly



Uro-Trapper Drain bag

- . Durabledrainbagsarenot made with natural rubber latex
- . Hasmeshfilterlargebagmouth, rapid drainage of fluids
- . 10/CS

Main features

737 0035 (Sterile)

5737.0036 (Non-sterile))

Size: 1630x20x510 mm



IG Tab-5737.0032

Uro-Trapper Support Ring

- It is usedwithdrainbagtoensure
- sanitation of surgery
- A flexiblesystem that attaches to any generalsurgical table tocontain fluids for hysteroscopy or cystoscopy procedures

Main features

can be installed on any kind of operating to

Soft, flexible band shape



Mattre



Maquet 1130 Split Leg

#508-0241 (8cm thickness) #508-0242 (10cm thickness)



Maquet 1130 Orthopedic

#508-0265 (8cm thickness) #508-0266 (10cm thickness)



MIZUHO 360

#508-0771 (8cm thickness) #508-0772 (10cm thickness)



MIZUHO 6700/6701

#508-0792 (8cm thickness) #508-0793 (10cm thickness)



Maquet Alphamaxx

#508-0289 (8cm thickness) #508-0290 (10cm thickness)



Maquet Alphastar 1132

#508-0313 (8cm thickness) #508-0314 (10cm thickness)



MIZUHO 3500

#508-0730 (8cm thickness) #508-0731 (10cm thickness)



MIZUHO 3501B

#508-0749 (8cm thickness) #508-0750 (10cm thickness)



MIZUHO5001/6001/6500/6600

#508-0692 (8cm thickness) #508-0693 (10cm thickness)



MIZUHO 3100

#508-0711 (8cm thickness) #508-0712 (10cm thickness)



Steris-Amsco 1080/2080

#508-0811 (8cm thickness) #508-0812 (10cm thickness)



Steris-Ams co1080/2080 w/Leg I

#508-0830 (8cm thickness) #508-0831 (10cm thickness)



Berchtold Operon D850/D750

#508-0198 (8cm thickness) #508-0199 (10cm thickness)



Maquet/Siemens 1130

#508-0219 (8cm thickness) #508-0220 (10cm thickness)



Steris-Amsco 3080/3085

#508-0849 (8cm厚度) #508-0850 (10cm厚度)



Steris-C-Max

#508-0868 (8cm thickness) #508-0869 (10cm thickness)









MANUFACTURING. DEVELOPMENT.CUSTOMIZATION

Qualität ist unsere Priorität

International group medical technology and electronics GmbH

Schongauer Str. 5 - 28219 Breme - Germany Tel: +49 (0) 421 8478 8804 www.ig-medical.de info@ig-medical.de



info@marinedasuministros.es (+34) 91 99 31 440 www.marinedasuministros.es

SEDE CENTRAL:

Rúa Agra da Bragua, nro. 6, bajo comercial. Barrio Ventorrillo 15010 A Coruña, Galicia, España.

DELEGACIÓN:

Calle de Velázquez, nro. 50, 5ta. planta. Venture X Velázquez 50 #1021 Barrio de Salamanca 28001 Madrid, España